


FRED PA-1

Simplicity as a priority



FRED PA-1

Dedicated to public access, SCHILLER's latest defibrillator, the FRED PA-1, has been designed so that even untrained users can save lives. The FRED PA-1 guides the rescuer step-bystep during the entire resuscitation process.

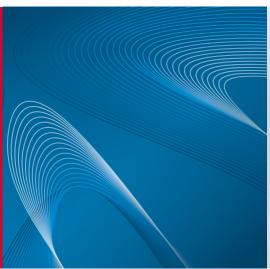





IT CONCERNS ALL OF US

More than 7 million people world-wide fall victim to a sudden cardiac arrest every year.

A sudden cardiac arrest can hit anyone, anywhere and at any time.

Defibrillation consists in delivering electrical energy to the heart in order to restore a proper heart rhythm. A defibrillator in combination with well-performed chest compressions is the only therapy that can save a patient's life.

The FRED PA-1 is a well-designed lifesaving equipment, which is highly visible and clearly recognizable in all environments.



HIGHEST USER-FRIENDLINESS

In emergency situations, procedures that would normally be clear and easy to follow can suddenly become confusing. Every detail counts and can make a difference in saving lives. That's why a defibrillator needs to be as simple as possible.

1-2-3 **STEPS**


Not only is the user interface streamlined with 1-2-3 steps, but spoken instructions and clear pictures also guide the user through the process.

NO DELAY


When you lift up the device cover, the FRED PA-1 starts up immediately. The electrodes are pre-connected: apply them to the patient's chest and let the FRED PA-1 guide you.

SHOCK DELIVERY POSSIBLE AS SOON AS A SHOCK IS RECOMMENDED

Resuscitation Guidelines emphasise the importance of minimizing the time-to-shock after Cardio Pulmonary Resuscitation (CPR). In order to reduce the hands-off time (time while no CPR can be performed), the FRED PA-1 features a fast analysis system that systematically pre-charges the correct energy so that the shock can be delivered as soon as it is recommended.

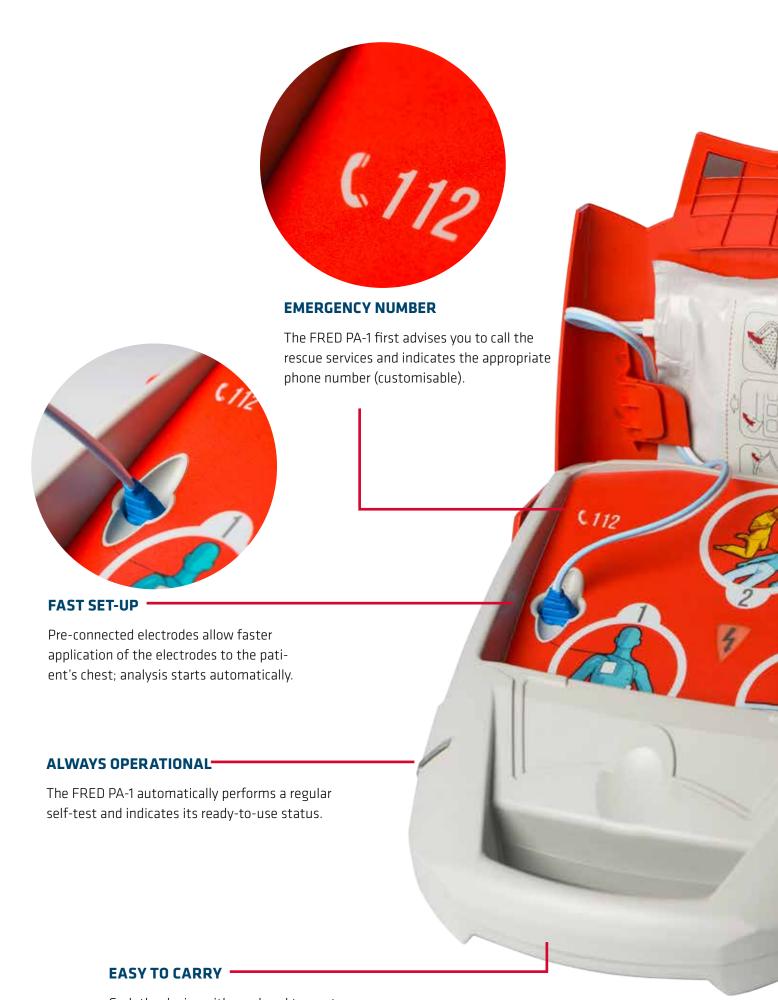

Grab the device with one hand to run to the scene of the incident.

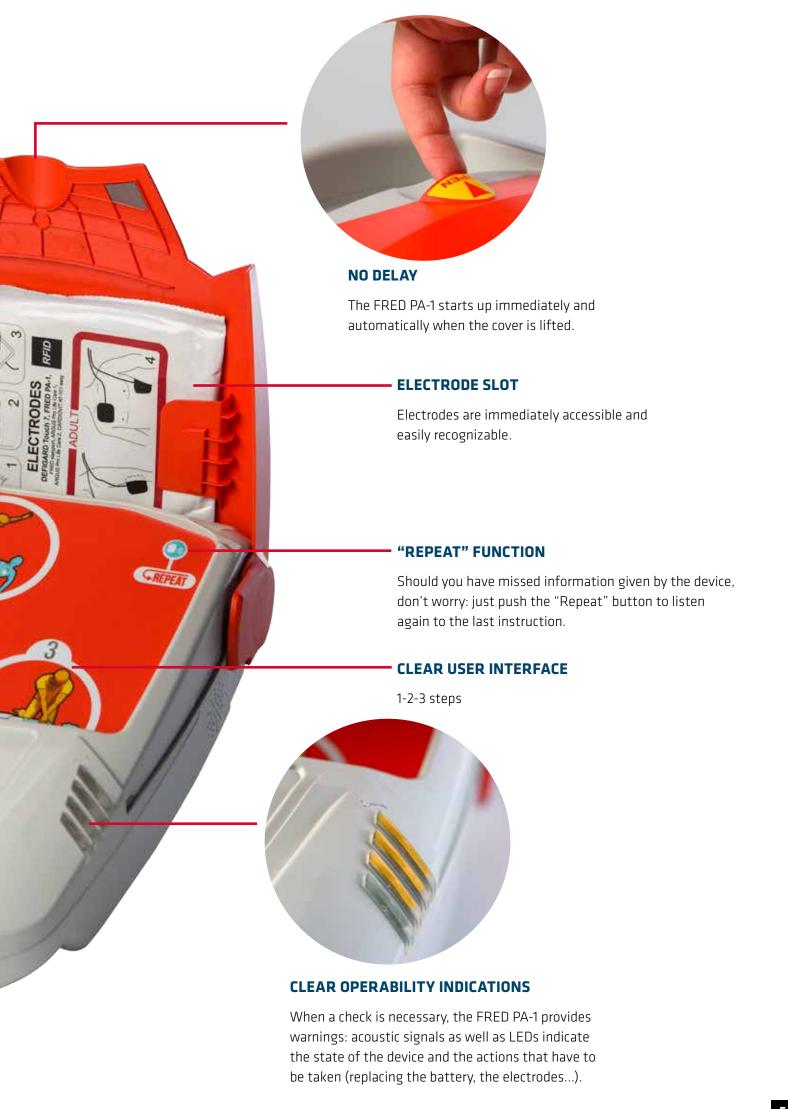



USED ANYWHERE

Sudden cardiac arrest can hit anyone, anywhere. The FRED PA-1 is designed for all locations, wherever you might need it.

With the ability to choose one of 3 languages¹, the FRED PA-1 can be installed in airports, hotels, railway stations, in any international environment.

With its high IP rating (IP55-protected against dust and water projection), the FRED PA-1 withstands even unfavourable conditions.

1 The three available languages can be customised at the time of order.

AUTOMATIC OR SEMI-AUTOMATIC

The FRED PA-1 analyses the patient's heart rhythm and determines whether a shock is necessary or not. Rescuers can fully rely on the result of the analysis, which is based on proven algorithms and decades of experience.

The FRED PA-1 is available either as semi-automatic or fully automatic device.

Delivering a shock to a victim of a sudden cardiac arrest should be done as soon as possible. In its automatic version, the FRED PA-1 delivers the shock without any action of the rescuer. In its semi-automatic version, the FRED PA-1 prompts the rescuer to deliver the shock by pressing the shock button.

DATA MANAGEMENT

The FRED PA-1 can store up to 2 hours of intervention data, including 500 events. Using a simple SD card, the retrieval of the intervention data as well as the software update are very easy to perform.

IT'S THE DETAILS THAT MAKE THE DIFFERENCE

EXTENDED WARRANTY PERIOD, UP TO 10 YEARS¹

SCHILLER provides peace of mind with up to 10 years of warranty – in the case of a warranty, the device is replaced with a new one.

READY WHEN YOU NEED IT

The FRED PA-1 has a long battery life (6 years) in order to be ready the day you need it. Battery level is checked during the self-test and a LED warns of low battery level.

ELECTRODE SHELF LIFE CHECK

Thanks to a highly innovative feature, the device is able to determine when the shelf life has expired and prompts the user to replace the electrodes.

1 Contact your local sales representative for more information

FOR ALL TYPES OF PATIENTS

Used with adult or children electrodes, the FRED PA-1 automatically recognises the electrode type and adjusts the defibrillation energy accordingly.

At the back of the device, a set of children electrodes can be safely stored.

CPR FEEDBACK

Defibrillation is not enough. In order to keep bringing oxygen to the organs, the blood flow must be kept. This is the purpose of chest compressions, and the reason why it's important to do it right. The FRED PA-1 helps you to achieve the highest efficiency in chest compressions by:

- indicating the right positioning of the hands
- setting a regular rhythm (metronome)
- providing a live feedback on the frequency of the chest compressions based on impedancemetry, without any additional sensor (FreeCPR).







SCHILLER Americas Inc. Doral, Florida 33172

North America:

Phone +1 786 845 06 20 Fax +1 786 845 06 02 sales@schilleramerica.com www.schilleramerica.com

Latin America & Caribbean

Phone +1 305 591 11 21 Fax +1 786 845 06 02 sales@schillerla.ch www.schillerla.ch



SCHILLER Asia-Pacific / Malaysia

52200 Kuala Lumpur, Malaysia Phone +603 6272 3033 Fax +603 6272 2030 sales@schiller.com.my www.schiller-asia.com



SCHILLER Handelsgesellschaft m.b.H.

A-4040 Linz Phone +43 732 70 99 Fax +43 732 757 000 sales@schiller.at www.schiller.at



China

Alfred Schiller (Beijing) Medical Equipment

100015 Beijing, China Phone +86 010 52007020 Fax +86 010 52007020-8016 info@schillerchina.com www.schillermedical.cn



Schiller medicinski instrumenti d.o.o.

10000 Zagreb Phone +385 1 309 66 59 +385 1 309 66 60 info@schillerzg.hr www.schiller.hr



France SCHILLER Médical S.A.S.

F-67162 Wissembourg Phone +33 3 88 63 36 00 +33 3 88 63 36 49 info@schiller.fr www.schiller-medical.com



France (distribution France) SCHILLER France S.A.S.

F-77608 Bussy St Georges Phone +33 1 64 66 50 00 +33 1 64 66 50 10 contact@schillerfrance.fr www.schiller-france.com



Germany SCHILLER Medizintechnik GmbH D-85622 Feldkirchen b. München Phone +49 89 62 99 81 0

Fax +49 89 609 50 90 info@schillermed.de www.schillermed.de



Hungary SCHILLER Diamed Ltd.

H-1141 Budapest Phone +36 1 383 4780 / +36 1 460 9491 Fax +36 1 383 4778 info@schillerhungary.hu www.schillerhungary.hu



SCHILLER Healthcare India Pvt. Ltd.

Mumbai – 400 059, India Phone +9122 6152 3333 / 2920 9141 +91 22 2920 9142 sales@schillerindia.com www.schillerindia.com



Poland

SCHILLER Poland Sp. z o.o.

PL-02-729 Warszawa Phone +48 22 843 20 89 / +48 22 647 35 90 +48 22 843 20 89 schiller@schiller.pl www.schiller.pl



www.schiller-cis.com

AO SCHILLER.RU 125124 Moscow, Russia Phone +7 (495) 970 11 33 Fax +7 (495) 956 29 10 mail@schiller.ru www.schiller.ru



Serbia SCHILLER d.o.o.

11010 Beograd Phone +381 11 39 79 508 +381 11 39 79 518 info@schiller.rs www.schiller.rs



Slovenia

SCHILLER d.o.o.

2310 Slovenska Bistrica Phone +386 2 843 00 56 Fax +386 2 843 00 57 info@schiller.si www.schiller.si



SCHILLER ESPAÑA, S.A. ES-28232 Las Rozas/Madrid Phone +34 91 713 01 76 +34 91 355 79 33 schiller@schiller.es www.schiller.es



Switzerland

SCHILLER-Reomed AG CH-8953 Dietikon Phone +41 44 744 30 00 Fax +41 44 740 37 10 sales@schiller-reomed.ch www.schiller-reomed.ch



Turkey SCHILLER TÜRKIYE Okmeydani-Sisli - Istanbul

Phone +90 212 210 8681 (pbx) +90 212 210 8684 info@schiller.com.tr www.schiller-turkiye.com



United Kingdom

SCHILLER ÜK Bellshill, ML4 3PR Phone +44 1698 744 505 +44 1698 744 474 sales@schilleruk.com www.schilleruk.com





